

CERAGEM

MASSAGE CHAIR

USER MANUAL

PAUSE M2
CGMEMCG-2101



CERAGEM[®]



Thank you for purchasing the CERAGEM Massage Chair.
Please read this manual carefully before starting any use.
Please store the manual where you can always find it after reading it.

This User's Manual contains the product warranty information.
This product is a home massage chair for household use. Do not use for commercial or industrial purposes.
We may improve the design and specifications of this product without any further notice.

This product is exclusively for 110V.
Please make sure to plug the power plug into a dedicated 110V outlet.
Do NOT use a 220V transceiver to 110V as it may result in safety hazards.



Users height whose exceed 6'3" or weight more than 220 lbs is not recommended to use the product.
Product may fall over and get damaged which may cause accidents or injuries to user.

CERAGEM
LIFESTYLE
MASSAGE CHAIR
USER MANUAL

PAUSE M2 USER MANUAL

MASSAGE CHAIR
CGMEMCG-2101

IMPORTANT SAFETY INSTRUCTIONS

When using an electrical appliance, basic precautions should always be followed, including the following:

Read all instructions before using this appliance:

DANGER - To reduce the risk of electric shock:

- 1) Always unplug this appliance from the electrical outlet immediately after using and before cleaning.

WARNING - To reduce the risk of burns, fire, electric shock, or injury to persons:


- 1) An appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
- 2) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- 3) Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
- 4) Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- 5) Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination and repair.
- 6) Do not carry this appliance by supply cord or use cord as a handle.
- 7) Keep the cord away from heated surfaces.
- 8) Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like.
- 9) Never drop or insert any object into any opening.
- 10) Do not use outdoors.
- 11) Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- 12) To disconnect, turn all controls to the off position, then remove plug from outlet.
- 13) Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
- 14) Connect this appliance to a properly grounded outlet only. See Grounding Instructions.

SAVE THESE INSTRUCTIONS

This appliance is intended for household use only

CONTENTS

- 01** Safety Precautions
- 02** Assembly Method
- 03** Name of each part and the components
- 04** How to use the remote controller
- 05** Care and Storage
- 06** Troubleshooting
- 07** FAQ
- 08** Product Specification

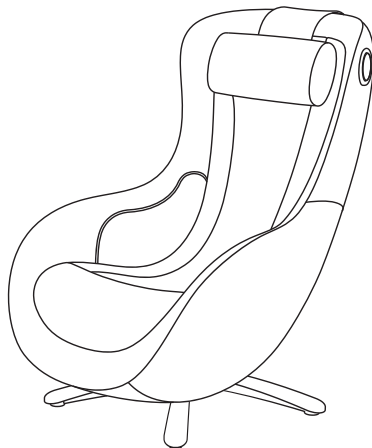

MASSAGE CHAIR
CGMEMCG-2101

01 Safety Precautions

The following safety precautions are intended to prevent unexpected danger or damage by ensuring safe and accurate use of the product.

When Installing

Remodeling, disassembly, or repairs should not be performed by anyone other than service personnel.

→ Can cause fire or electric shock.

Use the product on a flat surface.

→ Can cause product failure or an accident.

Regarding Power

If the product does not work or ceases to function, immediately disconnect the power plug from the outlet and contact the customer center for inspection and repair.

→ There is a risk of fire or electrocution due to electric shock or short circuiting.

If the cover of the product is torn or peeled off, it may cause a serious accident as clothes or hair can be caught inside. Ensure that there are no issues with the appearance of the product. If there is any damage, please contact the support center.

→ It may cause electric shock or injury.

Instruct children not to touch the power plug.

→ There is a risk of electric shock or injury.

Do not bend the power cord too hard and do not place heavy objects on it.

→ There is a risk of fire and electric shock.

Do not use multiple devices simultaneously when using a multi outlet adapter.

→ There is a risk of fire due to overheating of the outlet.

When unplugging the plug, pull it out of the outlet by grasping the plug.

→ If the power cord is severed, there is a risk of fire and electric shock.

Do not touch the power plug with wet hands, and when you clean or repair the product, disconnect the power plug from the outlet.

→ There is a risk of electric shock.

This product is exclusively for 110V. Please make sure to plug the power plug into a dedicated 110V outlet.

→ There is a risk of fire due to overheating of the outlet.

Be sure to use a grounding outlet when using an extension cord.

→ Using a non grounding outlet may cause electric shock or injury.

When there is a power outage, immediately turn off the power switch and unplug the power plug from the outlet.

→ Can cause electric shock or fire.

While operating the product, avoid pinching the AC cord or remote controller line between the center connection support or the four-foot support.

→ A short circuit may occur due to damage to the covering, so be careful when using it.

When using the product

Due to the risk of safety accidents, users with the following conditions should refrain from using the product:

The following persons should not use this product:

- People with thrombosis, embolism, severe venous, acute venous, etc.
- People with suspected acute cervical sprain (whiplash)
- People with symptoms of a sprain or sudden muscle contraction
- Pregnant women in early pregnancy or immediately after childbirth
- Those who have an abnormality in their spine due to an accident or disease in the past, or those with bent spine (scoliosis)
- Those who have an acute disease of the intestines or severe heart disorder
- Those who use an implantable electronic medical device such as implantable pacemakers
- Those who have a severe seizure symptom or psychogenic disorder (perceptual disorder)

-
- Those whose temperature sensation is lost due to severe peripheral circulation disorders such as diabetes, etc.
 - Those who have a risk of osteomyelitis or rickets
 - *Rickets: Symptoms of shrinking muscles throughout the body and an inability to stretch arms and legs
 - Those who are unable move their body freely
 - Those with or undergoing treatment for symptoms of disc disorders
 - Those receiving treatment at a medical institution for any disorder in their body

The following persons must consult a doctor before use:

- Those who have a risk of side effects due to the use of the product, e.g. pregnant women
- Those whose body temperature is higher than normal due to acute inflammation, chills, etc.
- Those with symptoms of acute skin disease, malignant tumor, or other skin disease
- Those who are receiving medical treatment
- Those who feel any bodily abnormalities
- Those who have been instructed by doctor not to undergo massage
- Those whose muscles have weakened or become lean due to old age
- Those who suffer from back pain

Do not use or store the product in a humid place such as a bathroom.

- It may cause electric shock, fire, failure, or mold.

Do not allow children or pets to play near or on top of the product and prevent them from using it.

- It may cause malfunction or injury.

Use under supervision if the user has a cognitive impairment or is inexperienced with the product.

- It may cause malfunction or injury.

Do not wear something hard on your head such as hair accessories during use.

- It may cause injury.

The operation time should not exceed a total of 30 minutes per day and 15 minutes at one time. Continuous use targeting the same area should be 5 minutes or less.

- Prolonged use may strain your muscles or nerves.

If you experience any of the following side effects while using the product, stop using it:

- Abrasions, strains, bruises, stinging/burning, or cuts caused by the product
- Damage to muscles, tendons, nerves, spine, etc. due to excessive use and pressure from the product
- Skin diseases and burns such as fractures, dislocations, sprains, etc.
- Respiratory and digestive disorders including asphyxiation, etc.

Do not use the product if you have any objects inside your pockets when massaging your hips or thighs.

- There is a risk of an accident or injury.

Make sure to use the product after checking that the power and remote control cords are neatly organized.

- It may cause an accident or injury or product failure.

Do not disconnect the power plug from the wall outlet or turn the power switch to [OFF] during the massage operation.

- This may result in an injury.

Do not use near open flames such as stoves or while smoking. Also, do not use it while mounted on or placed on top of a heating apparatus.

- It may cause fire.

Do not use after drinking alcohol.

- This may result in an accident or injury.

When operating product, please use the product by one person per product at a time. (Do not allow more than 2 or more persons to sit on the product)

- It may cause malfunction or injury.

Please do not use it on your bare skin.

- Direct massage on bare skin may cause skin damage.

Be careful not to get your legs caught in the remote control cable.


- This may result in an injury.

When cleaning or servicing the product, disconnect the power plug from the outlet.


- It may cause an electric shock or injury.

02 Assembly Method


1. Safely remove all the parts enclosed in the box.
(Please make sure that there are no parts left behind).
2. Place it on a flat surface to avoid scratching.
3. Lay the product forward, and then insert the center connecting support into the groove on the bottom of the product.
* Caution: Fit the first side of the 'center connecting support' to the fixing device.
4. Insert the 'four-foot support' all the way inside the fixing device.
5. Plug the connector to the sockets and turn the fastern to the right to fix it firmly.

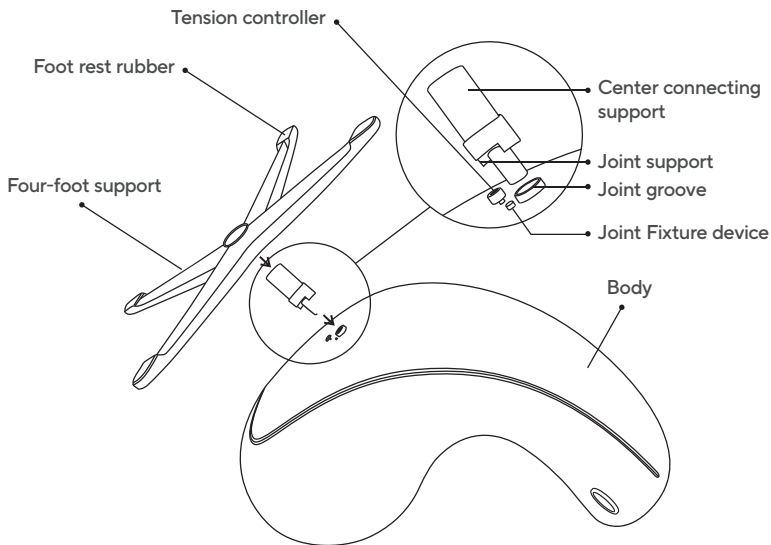


* When assembling and installing, put down the product slowly. If you rush this procedure, there is a risk of damage to the product or personal injury.

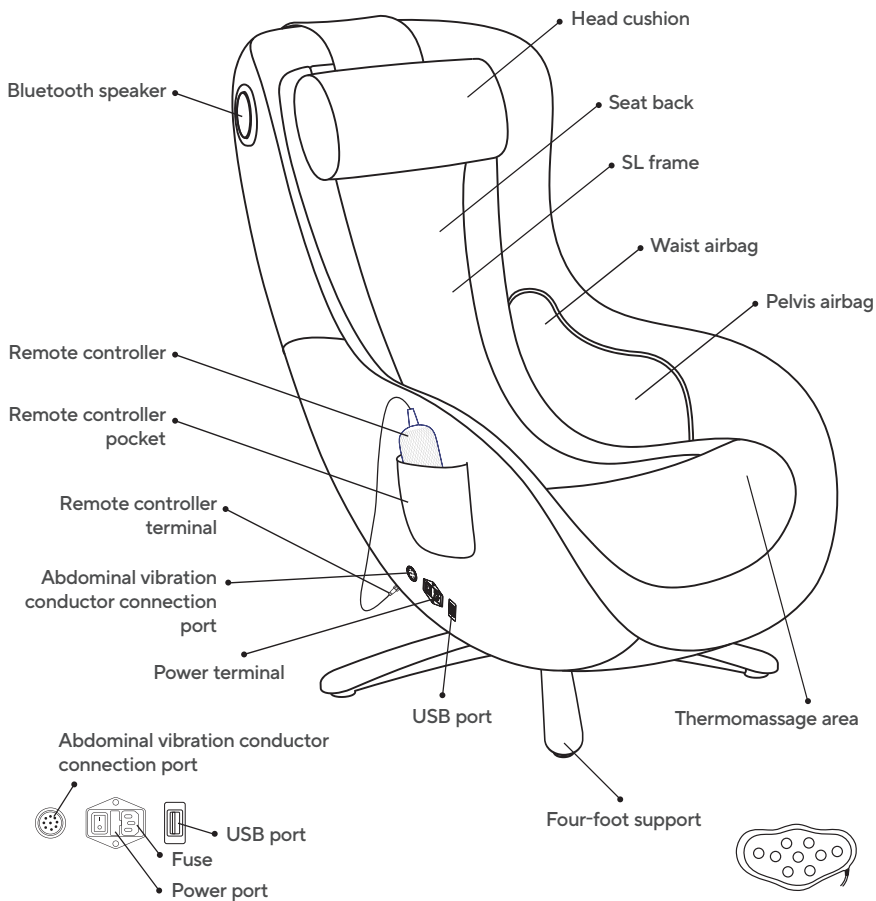

*Additional components: cushion,
1 electric abdominal ceramics

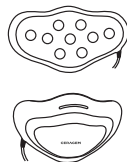

04 How to use the remote controller

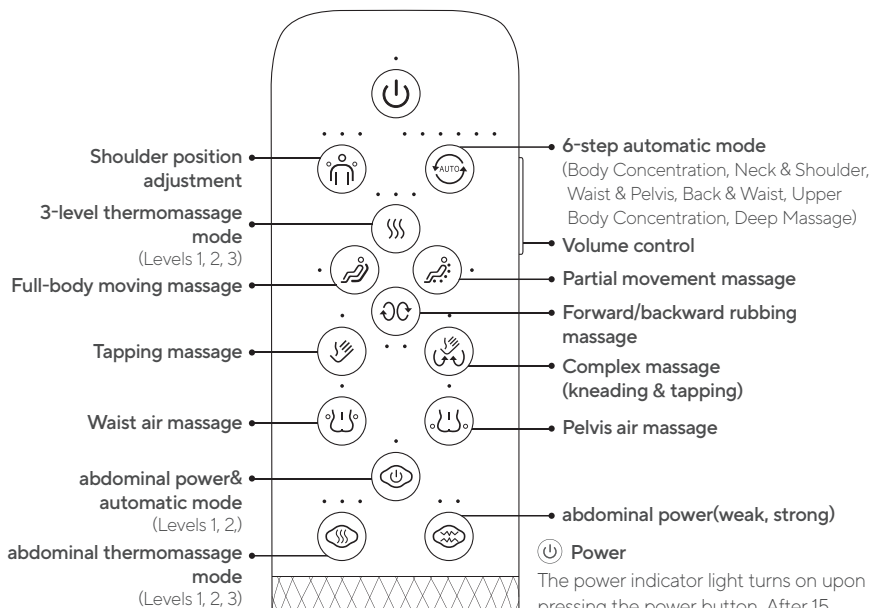

⏻ Power

The power indicator light turns on upon pressing the power button. After 15 minutes, the massage ends automatically and the massage ball returns to its position. After use the product, turn off the power and disconnect the power plug from the outlet to prevent malfunction.

🔄 6-step automatic mode

Each automatic mode lasts for 15 minutes, and after 'adjusting the shoulder position', press 'auto mode' to use it.

- 1 click: Full-body concentrating massage
- 2 click: Neck&Shoulder concentrating mode
- 3 click: Waist & Pelvic concentrating mode
- 4 click: Back & Waist mode
- 5 click: Upper body concentrating mode
- 6 click: deep massage mode

*Do not apply excessive force to the remote control or operate it with wet hands. If you operate the remote control with a pointed, hard object, the remote control may be damaged.

Range setting **Shoulder position adjustment button**

- 1 click: 5 cm higher
- 2 click: 5 cm lower
- 3 click: default position

 **Full body movement massage button**

- 1 click: Full-body moving massage
(Massage ball moves over the entire body)
- 2 click: Fixed massage
(Massage ball moves in its current position)

 **Partial movement massage button**

- 1 click: Partial movement massage
(The massage ball moves within a certain range)
- 2 click: Fixed massage
(Massage ball moves in its current position)

Abdominal massage function

 Press the power button to go directly to automatic mode.

- 1 click: vibration mode
- 2 click: Soft vibration mode

 **abdominal power (weak, strong)**

- 1 click: Fast
- 2 click: slow

 **abdominal thermomassage mode**

- 1 click: slow
- 2 click: medium
- 3 click: Fast

Function setting **Forward/backward rubbing massage**

- 1 click: Normal direction rubbing massage
- 2 click: Reverse direction rubbing massage
- 3 click: Stop

 **Tapping massage button**

- 1 click: Tapping
- 2 click: Fast tapping (1.5x)
- 3 click: Stop

 **Complex massage button**

- 1 click: Rubbing + Tapping
- 2 click: fast Rubbing + fast tapping (1.5 times)(1.5x)
- 3 click: Stop

Air massage button back **air massage**  **pelvis air massage**


- | | |
|-------------------|-------------------|
| 1 click: Lightly | 1 click: Lightly |
| 2 click: Strongly | 2 click: Strongly |
| 3 click: Stop | 3 click: Stop |

Bluetooth function

You can use the speaker by connecting a device with Bluetooth functionality to the main body

(ID: CGMEMCG-2101)

* If it is paired to another device, be sure to disconnect it before connecting.

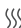 **3-level thermomassage mode**
Temperature control in 3 levels

- Level 1: slow
- Level 2: medium
- Level 3: Fast

When cleaning or maintaining the massage chair, be sure to disconnect the power plug from the outlet. Do not touch it with wet hands.

- There is a risk of an electric shock or injury.

Do not wipe with benzene thinner (diluent), alcohol or spray pesticides.

- It may cause discoloration, deterioration or fire.

Don't douse the product with water. Please wipe gently with a leather cleaner which is used for natural leather products.

Never iron the leather finish.

The leather finish is prone to color bleeding, so make sure to avoid wearing clothes susceptible to discoloration.

Continuous exposure to sunlight and ultraviolet rays may cause product discoloration. Avoid installation near windows or in areas with strong sunlight.

Information on natural leather sofas

PAUSE uses natural cow leather.

It is a natural phenomenon that each natural leather may vary in texture and color. Marks such as scars and wrinkles are rather characteristics of genuine leather.

Due to the features of the natural cowhide, the density of wrinkles and the texture of the skin may vary, so the surface pattern is not uniform, having slight variations.

- Leather is subject to different degrees of dyeing depending on its surface condition. The resulting slight color variation is natural appearance.

If any problem occurs while using the product, refer to the following solutions.

If the problem persists, turn off the power supply, disconnect the power plug from the outlet, and then ask the service center.

Users should never disassemble or repair the product at their own discretion.

CERAGEM is not responsible for any damage to the product caused by disassembly or repair at the user's discretion.

Problem	Solution
The product does not work.	<p>Check that the main power and power switch are turned on.</p> <p>Check the power plug connection state.</p> <p>Check that the correct voltage is used.</p>
You can hear a fricative sound from the massage machine part. You can hear the sound of the motor. Noise increases when using the tapping and kneading function.	Those noises are generated due to the mechanical structure and are not a malfunction. There may be noise even during normal operation.
During use, the remote control works, but the massage chair does not work.	<p>When the product is used for a long time, the massage chair automatically switches to overheating prevention mode due to increase in temperature.</p> <p>When the product is used for a long time, the massage chair automatically switches to overheating prevention mode due to increase in temperature. In this case, turn on the power and wait 30 to 50 minutes for the heat to cool.</p>
Both the remote control and the massage chair do not work.	<p>In case the working time has been set automatically, the set time has ended. Press the button again. Please turn the product off and on again.</p>

* If this does not work, please refer to the product warranty and contact customer support.

Q1. Can I get a massage on my back only?

After pressing the rubbing massage button, press the full-body moving massage button and press it again on the area where you want to be massaged. The massage ball will become fixed.

Q2. How many degrees does the chair lean back?

It is inclined up to 105-120 degrees when putting your feet on the ottoman, it will naturally lean back from the weight. (You cannot fix it while in the state of being inclined)

Q3. I can't connect to Bluetooth.

Turn the power off and then on. Then, delete CGMEMCG-2101 on the device to be connected, and reconnect it again.

*Bluetooth can only be connected to one device.

Q4. Will it create a lot of inter-floor noise?

Unlike other massage chairs, it barely creates any inter-floor noise because the main body is not attached to the floor.

In addition, you can further reduce the noise by placing a rug under the massage chair.

Q5. The operating sound seems louder than what I heard when trying it out at the store.

At the store, you may not be able to hear the sound of the massage chair due to the ambient noise. In a relatively quiet residential space, you may hear the machine sound when operating the massage chair, but this is a normal and not a cause for concern.

Q6. Can you transport the massage chair if I move to another place?

We do not provide a massage chair transport service.

The chair is not large in size and is lightweight, so it is easily moved and transported.

08 Product specification

Product name	PAUSE M2 (MASSAGE CHAIR)
Model name	CGM EMCG-2101
Size of the massage chair	W 25.6 in X D 44.1 in X H 39.7 in
Size of the ottoman	W 19.2 in X D 15.3 in X H 18.9 in
Net weight	81.3 lbs (massage chair), 9.2 lbs (ottoman)
Gross weight (including box)	About 112.4 lbs
Working voltage	100-240V, 50-60Hz
Power consumption	140W
Manufacture	Henglin Home Furnishings Co., Ltd.
Country of manufacture	China

Wait!

Warranty details can be found on www.ceragemus.com.



rtek
5000337



FCC ID 2ADIW-N-0002T

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Product name	PAUSE M2 (ABDOMINAL VIBRATION PROJECTOR)
Model name	CGM EMCG-2000
Size of the A.V.P	W 12.2 in X D 8.3 in X H 4.3 in
Net weight	3.1 lbs
Gross weight (including box)	About 3.6 lbs
Working voltage	DC 24V
Power consumption	48W
Manufacture	CERAGEM Co., Ltd.
Country of manufacture	Korea

Wait!

Warranty details can be found on www.ceragemus.com.



This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:


(1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

GROUNDING INSTRUCTION

This product must be grounded. If it should malfunction or breakdown, grounding provides a path of least resistance for electric current to reduce the risk of electric shock. This product is equipped with a cord having an equipment-grounding conductor and a grounding plug. The plug must be plugged into an appropriate outlet that is properly installed and grounded in accordance with all local codes and ordinances.

DANGER - Improper connection of the equipment-grounding conductor can result in a risk of electric shock. Check with a qualified electrician or serviceman if you are in doubt as to whether the product is properly grounded. Do not modify the plug provided with the product- if it will not fit the outlet, have a proper outlet installed by a qualified electrician.

This product is for use on a nominal 120-volt circuit, and has a grounding plug that looks like the plug illustrated in sketch A in following figure. A temporary adapter that looks like the adapter illustrated in sketch B and C may be used to connect this plug to a 2-pole receptacle as shown in sketch B if a properly grounded outlet is not available. The temporary adapter should be used only until a properly grounded outlet (sketch A) can be installed by a qualified electrician. The green colored rigid ear, lug, or the like extending from the adapter must be connected to a permanent ground such as a properly grounded outlet box cover. Whenever the adapter is used, it must be held in place by a metal screw.

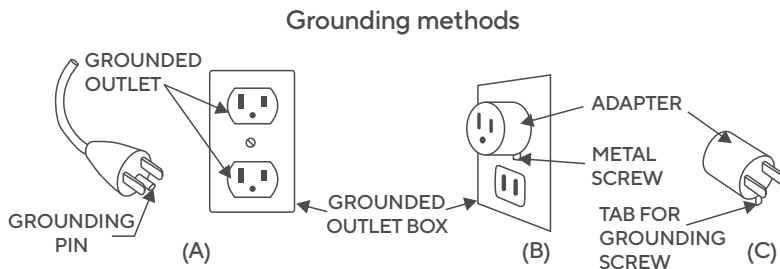

No adapter should be used with this product in Canada

CERAGEM®

Seller

CERAGEM CO., LTD.
10, Jeongja 1-gil,
Seonggeo-eup,
Seobuk-gu, Cheonan-si,
Chungcheongnam-do,
31045, Korea

Manufacturer

Henglin Home
Furnishings Co., Ltd.
No. 378, 380 Jiayi
Road, Dipu Street, ANJI
COUNTY Zhejiang
Province 313300, China

Customer Support Center

CERAGEM
INTERNATIONAL INC.
LOS ANGELES, CA 90010
3699 WILSHIRE BLVD.
STE 930
213-480-7070

FCC Requirement

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) this device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation. changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Please notice that if the FCC identification number is not visible when the module is installed inside another device, then the outside of the device into which the module is installed must also display a label referring to the enclosed module. This exterior label can use wording such as the following: "Contains FCC ID:2ADIW-N-0002T" any similar wording that expresses the same meaning may be used.

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with a minimum distance of 20cm between the radiator & your body. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

